



Gentaur Molecular Products BVBA
Address: Voortstraat 49, 1910 Kampenhout, Belgium
T: 0032 16 58 90 45 | **E:** info@gentaur.com
Websites: www.gentaur.com | www.maxanim.com

LSD

ELISA KIT INSTRUCTIONS PRODUCT # 130219 & 130215 (RTU)
FORENSIC USE ONLY

INTENDED USE: For the determination of trace quantities of LSD and/or other metabolites in human urine, blood, oral fluid.

DESCRIPTION

This ELISA (Enzyme-Linked ImmunoSorbent Assay) test kit is a qualitative one-step kit designed for use as a screening device for the detection of drugs and/or their metabolites. The kit was designed for screening purposes and is intended for forensic use only. It is recommended that all suspect samples be confirmed by a quantitative method such as gas chromatography/mass spectrometry (GC/MS).

ASSAY PRINCIPLES

This test kit operates on the basis of competition between the drug or its metabolite in the sample and the drug-enzyme conjugate for a limited number of antibody binding sites. First, the sample or control is added to the microplate. Next, the drug-enzyme conjugate is added and the mixture is incubated at room temperature. During this incubation, the drug in the sample or the drug-enzyme conjugate binds to antibody immobilized in the microplate wells. After incubation, the plate is washed to remove any unbound sample or drug-enzyme conjugate. The presence of bound drug-enzyme conjugate is recognized by the addition of K-Blue® Substrate (TMB). After a 30 minute substrate incubation, the reaction is halted with the addition of an acid stop. The test can be read visually or with a microplate reader equipped with a 450 nm filter. The extent of color development is inversely proportional to the amount of drug in the sample or control. In other words, the absence of the drug in the sample will result in a dark yellow color, whereas the presence of the drug will result in light yellow to no color development.

STORAGE AND STABILITY

This kit can be used until the expiration date on the label when stored refrigerated at 2-8°C. Store controls frozen if not used within 10 days. Note: Some kits require controls to be stored frozen immediately upon receipt. Reference kit label for details.

MATERIALS PROVIDED – SINGLE KIT (96 WELL)

1. **EIA Buffer:** 30 mL (ready-to-use). Phosphate buffered saline solution with bovine serum and a preservative. Provided for dilution of samples.
2. **Wash Buffer Concentrate (10X):** 20 mL. Phosphate buffered saline solution with a surfactant. Dilute 10 fold with deionized or ultrapure water before use. Diluted wash buffer is used to wash all unbound conjugate and samples from the plate after the conjugate incubation.
3. **K-Blue Substrate:** 20 mL (ready-to-use). Stabilized 3, 3', 5, 5' Tetramethylbenzidine (TMB) plus Hydrogen Peroxide (H₂O₂) in a single bottle. It is used to develop the color in the wells after washing. Light sensitive.
4. **Drug-Enzyme Conjugate:** 14 mL (ready-to-use). Drug-horseradish peroxidase conjugate. Do not dilute.
5. **Antibody Coated Plate:** A 96 well Costar plate, in strips of 8 break-away wells, coated with anti-drug antibody. The plate is ready for use as is. Do not wash.
6. **Acid Stop Solution:** 14 mL (ready to use). 1 N H₂SO₄ used to stop the enzyme reaction.
7. **Qualitative QC Positive Control:** 750 µL provided (synthetic human urine). Do not dilute.
8. **Qualitative QC Negative Control:** 750 µL provided (synthetic human urine). Do not dilute.

MATERIALS PROVIDED – BULK KIT (480 WELL)

1. **EIA Buffer:** 200 mL (ready-to-use). Phosphate buffered saline solution with bovine serum and a preservative. Provided for dilution of samples.
2. **Wash Buffer Concentrate (10X):** 100 mL. Phosphate buffered saline solution with a surfactant. Dilute 10 fold with deionized or ultrapure water before use. Diluted wash buffer is used to wash all unbound conjugate and samples from the plate after the conjugate incubation.
3. **K-Blue Substrate:** 100 mL (ready-to-use). Stabilized 3, 3', 5, 5' Tetramethylbenzidine (TMB) plus Hydrogen Peroxide (H_2O_2) in a single bottle. It is used to develop the color in the wells after washing. Light sensitive.
4. **Drug-Enzyme Conjugate** (ready-to-use): 5 X 14 mL. Drug-horseradish peroxidase conjugate. Do not dilute.
5. **Antibody Coated Plate:** 5 X 96 well Costar plates, in strips of 8 break-away wells, coated with anti-drug antibody. The plates are ready for use as is. Do not wash.
6. **Acid Stop Solution:** 90 mL (ready-to-use). 1 N H_2SO_4 used to stop the enzyme reaction.
7. **Qualitative QC Positive Control:** 5 X 750 μ L (synthetic human urine). Do not dilute.
8. **Qualitative QC Negative Control:** 5 X 750 μ L (synthetic human urine). Do not dilute.

MATERIALS NEEDED BUT NOT PROVIDED

1. Deionized water.
2. Precision pipettes that range from 10 μ L - 1000 μ L and disposable tips.
3. Graduated cylinder to dilute and mix wash buffer.
4. Plate cover or plastic film to cover plate during incubation.
5. Microplate reader with a 450 nm filter.
6. Cut-off calibrator.

OPTIONAL MATERIALS

1. Microplate shaker.

PRECAUTIONS AND NOTES

1. **DO NOT** use kits or components beyond expiration date.
2. **DO NOT** mix conjugates and plates from different kit lots.
3. **DO NOT** pipette reagents by mouth.
4. Pour K-Blue Substrate out of the bottle into a clean reservoir. To prevent contamination of the substrate, **DO NOT** pipette out of the bottle.
5. All specimens should be considered potentially infectious. Exercise proper handling precautions.
6. Keep plate covered except when adding reagents, washing or reading.
7. Kit components should be refrigerated at all times when not in use.
8. Keep the controls frozen if storing longer than 10 days. Avoid repeated freeze-thaw cycles. Note: Some kits require controls to be stored frozen immediately upon receipt. Reference kit label for details.
9. Use aseptic technique when opening and removing reagents from vials and bottles.
10. **DO NOT** smoke, eat or drink in areas where specimens or reagents are being handled.
11. Do not substitute DI water for the wash step of this protocol. Use only our wash buffer.
12. Do not use Sodium Azide with samples, standards and/or calibrators. Concentrations as low as 0.01% Sodium Azide will reduce the color of this assay.
13. Do not reuse wells, they are for one use only.

PROCEDURAL NOTES

1. Desiccant bag must remain in foil pouch with unused strips. Keep ziplock pouch sealed when not in use to maintain a dry environment.
2. Use clean pipette tips for the buffer, drug-enzyme conjugate, controls and samples.
3. Before pipetting a reagent, rinse the pipette tip three times with that reagent.
4. When pipetting into the wells, **DO NOT** allow the pipette tip to touch the inside of the well or any of the reagent already inside the well. This may result in cross contamination.
5. Controls and samples should be assayed in duplicate.
6. Before substrate addition, wipe the outside bottom of the wells with a lint-free wiper to remove dust and fingerprints.
7. Gently mix specimens and reagents before use. Avoid vigorous agitation.

SAMPLE TREATMENT

Recommended minimum sample dilutions are listed below. These dilutions may change based on your laboratory's determination. All sample dilutions should be made in the provided EIA Buffer.


- a. **Urine:** Samples should be assayed neat. This assay is optimized at a 0.5 ng/mL cutoff when human urine samples are assayed neat. Samples containing sodium azide at a concentration of 0.01% or higher should be avoided unless laboratory calibrators contain the same level of sodium azide.
- b. **Whole blood:** Samples should be diluted 5X with EIA Buffer. When diluted 5X, cutoffs ranging from 2-5 ng/mL may be used.
- c. **Oral fluid:** Samples may be assayed neat.

TEST PROCEDURES

The following test procedures can be run manually or on an automated instrument.

1. Determine the number of wells to be used.
2. Gently mix the ready to use conjugate solution by inversion. Do not vortex.
3. Add 50 μ L of sample, laboratory calibrators and controls to the appropriate wells in duplicate. DO NOT dilute positive and negative controls.
4. Add 100 μ L of the ready to use drug-enzyme conjugate to each well. Use 8-channel pipetter or 12-channel pipetter for rapid addition.
5. For manual runs, mix by gently shaking plate. A microplate shaker may be used.
6. Cover plate with plastic film or plate cover and incubate at room temperature for 45 minutes.
7. During the conjugate incubation, dilute concentrated wash buffer 10 fold with deionized water (i.e. 20 mL of concentrated wash buffer plus 180 mL of deionized water). Mix thoroughly. Diluted wash buffer is stable for 5 days at room temperature or 7 days at 2-8°C.
8. Once the incubation is complete, dump or aspirate the liquid from the wells. Tap the plate on a clean lint-free towel to remove any remaining liquid in the wells.
9. Wash each well with 300 μ L of diluted wash buffer. Manual Wash: For manual wash procedures repeat for a total of 3 washings, invert and tap dry the plate following each step. After completing the last wash step wipe the bottom of the wells with a lint-free towel to remove any liquid on the outside of the wells. Automated Wash: If an automated plate washer is used wash the plate for a total of 5 washings with 300 μ L of diluted wash buffer. It is important for the automated washer to conduct a final aspirate cycle to eliminate residual amounts of wash buffer. Residual amounts of buffer in the wells will affect assay performance. Note: DI water should never be used for the plate wash.
10. Add 100 μ L of the K-Blue Substrate to each well. For manual runs, use a multi-channel pipetter for best results.
11. Incubate at room temperature for 30 minutes.
12. Add 100 μ L of the Acid Stop (1N H₂SO₄) to each well to stop enzyme reaction. Mix gently before measuring absorbance. For automated systems a 10 second shake is sufficient. Measure the absorbance at a wavelength of 450 nm. Wells should be read within 2 hours of stopping the reaction.

STANDARD CURVE IN EIA BUFFER

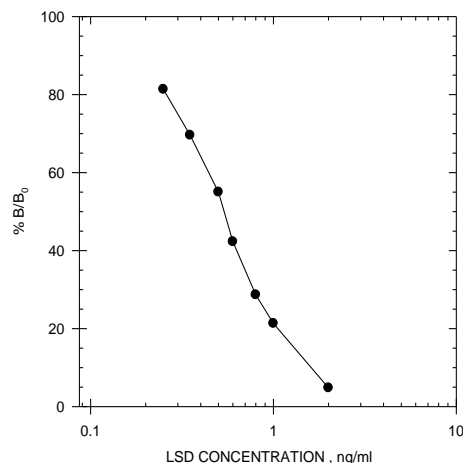

SENSITIVITY

Compound	I-50 in EIA Buffer (ng/mL)	I-50 in Human Urine (ng/mL)
LSD	0.63	0.42

The term I-50 is used to define the sensitivity of the test. This number is derived from a standard curve generated with the drug in EIA Buffer. The drug concentration that shows 50% less color activity than the zero standard is considered to be the I-50.

SPECIFICITY

Compound	Compound Concentration (ng/mL)	LSD Equivalents (ng/mL)	% Cross- Reactivity
LSD	0.54	0.54	100
LAMPA	0.35	0.54	153
Fentanyl	0.97	0.54	56
Nor-/Nor-iso-LSD	1.48	0.54	36
Acetyl Fentanyl	4.15	0.54	13
Iso-LSD	21.1	0.54	2.6
2-oxo-3-OH LSD	26.2	0.54	2.1
2-oxo-3-OH LAMPA	45.5	0.54	1.2
Norfentanyl	56.8	0.54	1.0
Lobeline	675	0.54	0.08
Clenbuterol	769.7	0.54	0.07
Methylene Blue	948.7	0.54	0.06
Procaineamide	2175.2	0.54	0.02

Note: LSD equivalents represents 50% B/B₀ assay displacement in EIA Buffer.

The compounds having cross-reactivity below 0.01% did not show any significant reaction up to 10µg/mL.

ALL THE FOLLOWING HAVE A CROSS-REACTIVITY <0.01%.

Acepromazine; Acetaminophen; Acetylsalicylic Acid; E-amino-n-caproic Acid; Amitriptyline; Ascorbic Acid; Benzoic Acid; Caffeine; Chlorodiazepoxide; Chlorpromazine; Codeine; Cotinine; Dexamethasone; Dextromethorphan; Diclofenac; Dimethyl Sulfoxide; Dipyrone; Doxepin; Ephedrine; Erythromycin; Ethyl p-amino-benzoate; Fenoprofen; Flunixin; Folic Acid; Folinic Acid; Furosemide; Gemfibrozil; Gentsic Acid; Glipizide; L-Glutamic Acid; Glutethamide; Glycopyrrolate; Heparin; Hippuric Acid; Hordenine; Hydrocortisone; Ibuprofen; Imipramine; Isoxsuprine; JWH-018; Lidocaine; Meperidine; Mescaline; Metaproterenol; Methadone; Methaqualone; Methocarbamol; Methylprednisolone, Nalorphine; Naproxen; Niacinamide; Nicotine; Nortriptyline; Orphenadrine; Oxyprenbutazone; Penicillin G-Potassium; Penicillin G-Procaine; Pentoxifylline; Phenothiazine; Phenylbutazone; Polyethylene Glycol; Prednisolone; Primadone; Procaine; Promazine; Pseudoephedrine; Pyrantel; Pyrimethamine; Pyrilamine; PCP; Quinidine; Quinine; Salbutamol; Salicylamide; Salicylic Acid; Salvinorin A; Theophylline; Thiamine; Trimethoprim; Trimipramine; Uric Acid

RESULTS INTERPRETATION

Positive Result: Samples with an absorbance less than or equal to the laboratory's designated cutoff calibrator should be considered positive. All positive samples should be confirmed by a quantitative method such as GC/MS.

Negative Result: Samples with an absorbance greater than the laboratory's designated cutoff calibrator should be considered negative.

Note: The kit was designed for screening purposes only. It is recommended that all suspect samples be confirmed by a quantitative method such as GC/MS or HPLC.

Products labeled **FORENSIC USE ONLY** are tests that have not been FDA approved.